

**INDUCTIVELY COUPLED
PLASMA EMISSION
SPECTROSCOPY (ICP-
OES) - TORONICP™
700TP
(PETROCHEMICAL
EDITION)**



Inductively Coupled Plasma Emission Spectroscopy (Petrochemical Edition) is a specialized inductively coupled plasma spectrophotometer engineered for high-accuracy elemental analysis in demanding petrochemical environments. Using advanced inductively coupled plasma emission spectroscopy, the system generates a stable, high-temperature inductively coupled plasma that enables reliable, simultaneous multi-element detection in fuels, lubricants, additives, and organic chemical products.

This petrochemical-focused platform incorporates optimized sample introduction, oxygen assisted combustion, and low-temperature fog chamber technology to effectively process complex organic matrices with minimal preparation. Intelligent automation and robust plasma control allow consistent, laboratory-grade results even for users without prior experience in inductively coupled plasma analysis, making the system well suited for routine quality control and advanced petrochemical testing laboratories.



APPLICATIONS

Applications

The Inductively Coupled Plasma Emission Spectroscopy (Petrochemical Edition) ToronICP™ 700TP is purpose-built for petrochemical and energy-related elemental analysis, including:

- Fuel analysis of gasoline and diesel with direct injection, eliminating complex digestion steps
- Lubricating oil and additive testing, measuring over 20 wear metals and additive elements
- Used oil analysis for condition monitoring and contamination assessment
- Petrochemical raw material testing, including methanol and ethanol
- Automotive fluid analysis, such as urea solutions and specialty oils
- Oilfield and tracer analysis, including rare earth element determination





FEATURES

Inductively Coupled Plasma Emission Spectroscopy (ICP-OES) - ToronICP™ 700TP (Petrochemical Edition) Key Features

- Inductively Coupled Plasma Emission Spectroscopy (ICP-OES) for fast, simultaneous multi-element analysis
- STCC-BOX Refrigerated Fog Chamber, providing rapid cooling to $-25\text{ }^{\circ}\text{C}$ for stable volatile sample introduction
- Direct organic sample injection, reducing preparation time and improving workflow efficiency
- Oxygen-assisted injection technology, minimizing carbon deposition and improving plasma stability

- High-sensitivity atomization system, optimized for petrochemical matrices
- Precision temperature control, enhancing analyte stability and repeatability
- User-friendly automated operation, suitable for both expert laboratories and first-time ICP users

THEORY & METHOD

Theory and Working Principle

Inductively Coupled Plasma Optical Emission Spectroscopy (ICP-OES) is a powerful method for elemental analysis. This technique utilizes an extremely hot, stable plasma, generated by energizing argon gas with a radio-frequency electromagnetic field, which can reach temperatures exceeding 10,000 K.


The process involves converting liquid samples into a fine aerosol and introducing them into this high-temperature plasma. Within the plasma, the elements in the sample become excited and subsequently emit light at specific wavelengths unique to each element. An inductively coupled plasma spectrophotometer then detects this emitted light, allowing for the simultaneous identification and precise quantification of multiple elements.

ICP-OES is highly valued for its speed, accuracy, high sensitivity, minimal chemical interference, and broad linear measurement range. These characteristics make it a reliable and versatile tool across various fields, including environmental monitoring, metals analysis, industrial quality control, and scientific research.

TECHNICAL SPECIFICATIONS

Inductively Coupled Plasma Emission Spectroscopy (Petrochemical Edition) ToronICP™ - 700TP Technical Specifications

Model	ToronICP™ -700TP
Analytical Technique	Inductively Coupled Plasma Emission Spectroscopy (ICP-OES)
Application Focus	Petrochemical, fuels, lubricants, additives, organic chemicals
Sample Introduction	Direct injection with high-efficiency atomizer
Fog Chamber	STCC-BOX refrigerated fog chamber
Fog Chamber Temperature	Down to –25 °C
Cooling Time	Stable temperature achieved in ~10 minutes
Plasma Support	Oxygen-assisted injection system
Elements Analyzed	Metals and trace elements in organic and oil-based samples
Operation Mode	Fully automated
Sample Preparation	Minimal or no digestion required
Intended Use	Laboratory and industrial petrochemical analysis
Analytical speed	10 elements every minute
Wavelength range	195-800nm
Accuracy	Relative standard deviation RSD≤2%
Single package size	1680 x 740 x 1320 mm
Weight	100 KG

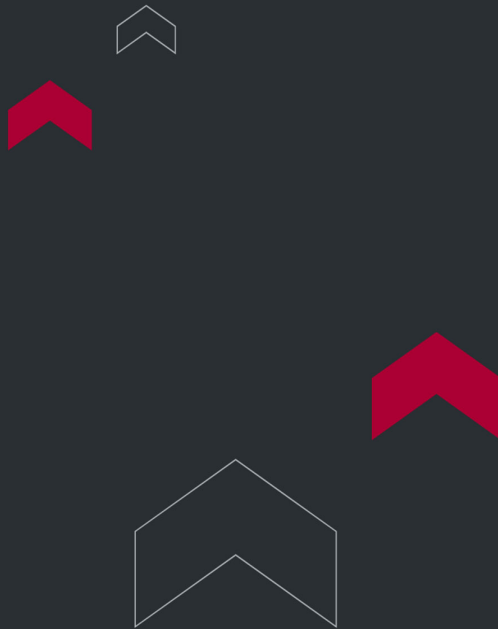

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

