



HOSPITAL ICU/CCU BED - TORONCARE 1070

The ToronCare™ 1070 is an advanced ICU/CCU electric hospital bed engineered for critical care units, post-surgical recovery, and intensive patient monitoring. Built with durable steel framework and premium ABS & PP plastic, this adjustable medical bed features a 10.2-inch touchscreen interface, integrated weight scale, and off-bed alarm system for enhanced patient safety and precise positioning with 250 kg load capacity.

This electric adjustable hospital bed offers user-friendly touchpad controls for back-rest (0-65°), leg-rest (0-40°), Trendelenburg (0-10°), and hi-low adjustment (530-710 mm). Advanced features include simultaneous back-rest and leg-rest adjustment, patient lock-out system, backup battery for emergency power failure, and reset button for quick flat positioning. The custom bed board enhances patient comfort, while electric central brake castors and detachable PP head/foot boards ensure safe operation. Optional X-ray radiolucent sections enable diagnostic imaging without patient transfer.

Ideal for rehabilitation and long-term critical care, the ToronCare™ 1070 delivers comprehensive patient monitoring capabilities and caregiver efficiency through intuitive touchscreen controls and data storage functionality for clinical environments.

APPLICATIONS

- **Critical and Intensive Care Units** – Provides advanced electric controls and touchscreen interface for precise patient positioning.
- **Post-Surgical Recovery** – Supports efficient adjustment of back-rest and leg-rest simultaneously for patient comfort.
- **Patient Monitoring and Safety** – Features weight scale, off-bed alarm, and lock-out system for enhanced security.
- **Emergency Preparedness** – Equipped with a backup battery and reset button to ensure quick response during power failure.
- **Diagnostic and Imaging Support** – Optional X-ray radiolucent sections allow imaging without patient transfer.
- **General Hospital Wards** – Durable, user-friendly design ensures long-term use and easy caregiver operation.
- **Rehabilitation and Long-Term Care** – Offers patient comfort with custom bed board, smooth adjustments, and ergonomic positioning.

FEATURES



- Built with durable steel framework and high-quality ABS & PP plastic
- Controlled by user-friendly touchpad (replacing traditional keypads, which have limited service life due to mechanical failure)
- 10.2 inch touchscreen at the foot board displays all the bed posture positions and features a patient lock-out system
- Custom bed board improves softness and bending of the bed, increasing the comfort of patients
- Electrically controlled bed positions provide easy adjustment of patient posture, with an indicator light to show the work status
- Back-rest and leg-rest can be adjusted up/down at the same time for more efficient operation
- Bed position and weight scale data can be stored for more efficient operation
- Off-bed alarming
- Reset button adjusts the bed to the flat / lowest setting from any position
- Backup battery provides protection for emergency power failure
- Detachable head and foot board made from high-quality PP plastic
- 4 pcs drop-down PP protective railing, controlled by gas spring
- Central brake castors, brake and un-brake controlled by electric operation
- IV pole mount bracket at all 4 corners of the bed frame

Standard Accessories (Included)

- IV Pole


Optional Accessories (Not Included)

- Monkey pole
- X-ray radiolucent head and back section

TECHNICAL SPECIFICATIONS

Dimensions

Overall dimensions (LxWxH)	2260mmx1050mm x 530-710mm 89.0in x 41.3in x 20.9-28.0in
Hi-Low	710mm-530mm 28.0-20.9in
Bed Board dimensions (LxW)	1940mmx900mm 76.4in x 35.4in
Diameter of central brake castors	150mm 5.9in
Bed Articulations	
Back-rest	0-65o
Leg-rest	0-40o
Trendelenburg and Rev. Trendelenburg	0-10o
Weight	
Net Weight	158 Kg 348 lbs
Load Capacity	<250 Kg <551 lbs

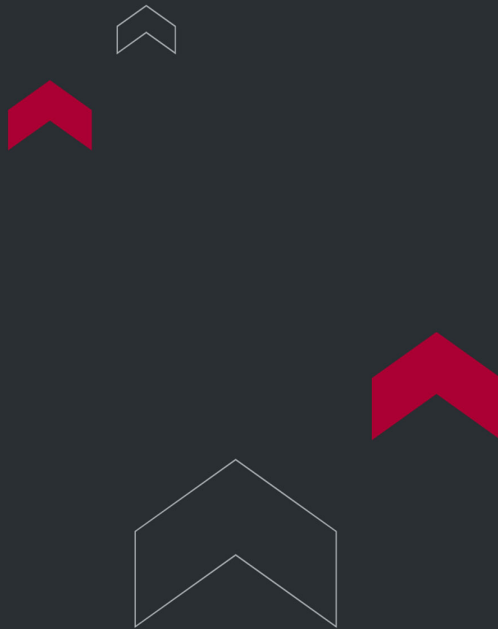

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 🇨🇦

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 🇺🇸

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 🇦🇪

Phone: +971 488 19252
Email: me@torontech.com

