



UV WEATHERING TEST CHAMBER

UV Weathering Test Chamber replicates the harmful effects of sunlight, rain, and dew on materials. It condenses months or years of natural weathering into a matter of days or weeks. The testing process alternates between UV light exposure and humidity cycles, with materials subjected to carefully regulated high-temperature conditions.

APPLICATIONS

Fluorescent UV lamps simulate sunlight, while condensation and water spray mimic dew and rain. This process ensures accurate and efficient weathering simulations for assessing material durability.

Real Simulation

The UV Weathering Test Chamber is equipped with UVA: 340nm lamps, which accurately replicate the short-wavelength spectrum of sunlight. This segment of sunlight is particularly damaging to materials like plastics and paints. This precise simulation ensures reliable results, closely correlating with outdoor exposure conditions.

Cost-Effective Design

Both purchase and operation costs are remarkably economical. The tower-style UV test chamber features an efficient design, using budget-friendly fluorescent UV lamps to mimic sunlight and standard tap water for condensation.

Effortless Operation

With a compact and straightforward design, the chamber is easy to install, operate, and maintain.

1. **User-friendly interface with intuitive programming**
2. **Continuous monitoring and display of exposure conditions**
3. **Built-in diagnostic alarms and maintenance alerts**

Extensive Expertise

Our advanced aging test chambers and services are crafted with extensive experience in standardized weathering testing. Our team actively contributes to ISO, ASTM, IEC, GB, and other global organizations, shaping best practices for material durability assessments.

UV Weathering Test Chamber

Wet Environment Simulation

The tower-style UV Test Chamber offers two moisture simulation methods: condensation and water spray. While condensation is the primary approach for replicating outdoor moisture damage, the chamber also supports a water spray cycle to meet specific application requirements. This dual functionality ensures a comprehensive simulation of real-world moisture conditions.

Condensation Cycle

In many outdoor settings, materials remain wet for up to 12 hours daily. Studies reveal that dew, rather than rain, is the leading cause of outdoor wetness. To replicate this, the tower UV Test Chamber uses a unique condensation feature. During this cycle:

- Water vapor is generated by heating a water tank at the chamber's base.
- Hot water vapor fills the chamber, creating 100% relative humidity at elevated temperatures.
- The test sample forms the chamber's side wall, with its back exposed to ambient air.

The cooling effect of the surrounding air reduces the sample surface temperature slightly below that of the water vapor. This temperature difference ensures continuous condensation of water vapor on the sample surface.

The condensed water is stable, pure distilled water. Even when ordinary tap water is used, the water collected on the sample surface remains free of impurities, delivering reliable testing conditions for material durability assessments.

Condensation: Reliable and Efficient Moisture Simulation

- **Natural Consistency:** Replicates real-world moisture conditions effectively.
- **Accelerated Testing:** Speeds up the process by raising temperatures.
- **Tap Water Usage:** Operates efficiently using standard tap water.
- **Pure Condensation:** Produces distilled water, ensuring no water spots or contamination on test samples.
- **Low Maintenance:** Designed for easy upkeep and hassle-free operation.

Water Spray: Simulating Thermal Shock and Mechanical Erosion

The water spray feature replicates real-world conditions where sudden rain showers cool sun-heated materials, causing rapid temperature changes. This process tests materials for their ability to withstand thermal shock. Additionally, the water spray simulates the mechanical erosion caused by rainfall, assessing the durability of materials against such wear and tear.

- **Thermal Shock Simulation:** Mimics abrupt temperature fluctuations from rain on heated surfaces.
- **Mechanical Erosion Testing:** Evaluates the impact of rainfall on material integrity.
- **Pure Water Requirement:** Unlike the condensation system, water spray requires laboratory-grade pure water for optimal performance.

UV and Sunlight Simulation

Ultraviolet light is a leading cause of photodegradation in outdoor products. The UV Weathering Test Chamber uses advanced fluorescent UV lamps to replicate the critical short-wave UV light found in sunlight, effectively simulating the aging

effects on materials caused by prolonged exposure to light.

UV Lamp Options

The chamber offers a range of high-quality UV fluorescent lamps to accommodate various testing conditions. These lamps feature:

- **Stable Spectrum Output:** Ensures consistent performance throughout their lifespan.
- **Rigorous Quality Control:** Each batch undergoes thorough testing to meet exacting standards before being approved for sale.

This combination of precise simulation and reliable components ensures accurate and repeatable testing results.

UV Weathering Test Chamber

UVA Lamps: Precision Sunlight Simulation

UVA lamps provide accurate replication of sunlight in critical UV regions, making them essential for material durability testing:

1. **UVA340 Lamps:** Simulate the short-wavelength UV spectrum of sunlight, specifically between 365nm and 295nm, effectively replicating the most damaging rays.
2. **UVA351 Lamps:** Designed to mimic UV light that passes through window glass, ideal for testing indoor exposure effects.
3. **Cool White Lamps:** Suitable for simulating artificial lighting conditions, such as those in office environments.

Each lamp undergoes a rigorous series of tests to ensure performance consistency and quality, delivering reliable results for a variety of testing scenarios.

UV Weathering Test Chamber

UVB Lamps: High-Intensity UV Testing Solutions

UVB lamps offer enhanced ultraviolet light for rigorous testing of material durability:

1. **UVB313EL Lamps:** Emit short-wave UV light stronger than natural sunlight at the Earth's surface, ideal for:
2. Quality control procedures.
3. Product development processes.
4. Testing materials with high weather resistance.
5. **UVB311EL+ Lamps:** Deliver exceptionally high irradiance levels, making them suitable for demanding fluorescent UV tests requiring intensified conditions.

These lamps provide reliable tools for simulating extreme UV exposure, ensuring precise evaluation of material performance.

UV Weathering Test ChamberUV Weathering Test Chamber

Irradiance Control: Ensuring Accurate and Consistent Testing

The intensity and spectrum of UV light directly influence the rate and type of material aging. Accurate control of irradiance is crucial for reliable and repeatable test results.

- **Closed-Loop Precision:** The irradiance control system features a precise closed-loop mechanism, automatically maintaining the desired light intensity using a feedback loop.
- **Continuous Monitoring:** The system continuously tracks UV light intensity and adjusts the lamp output to account for changes caused by factors like lamp aging.
- **Superior Accuracy:** Delivers better reproducibility and consistency compared to manual control systems.


This advanced control ensures that weathering tests provide consistent, reliable data for evaluating material performance.


UV Weathering Test Chamber

THEORY & METHOD

Real Simulation

Dew, not rain, is responsible for most of the humidity in outdoor exposures.

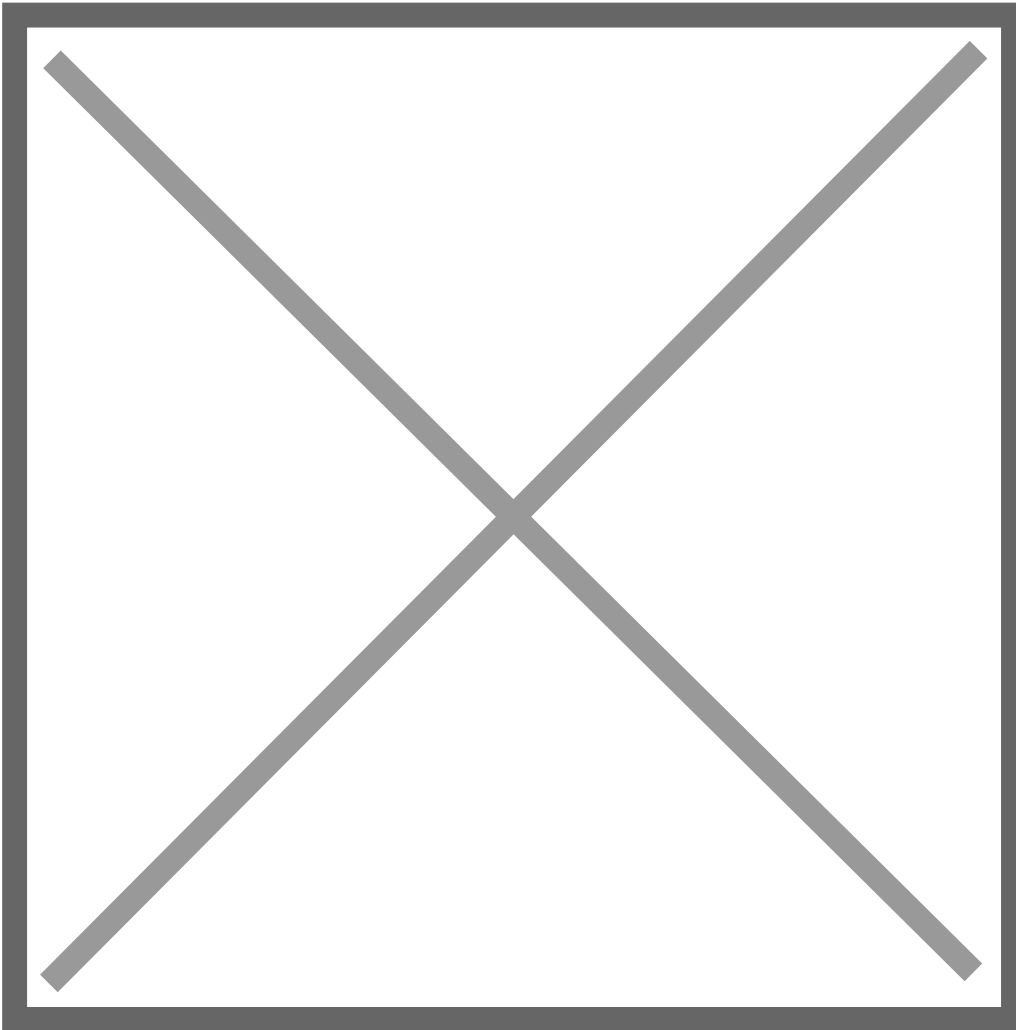

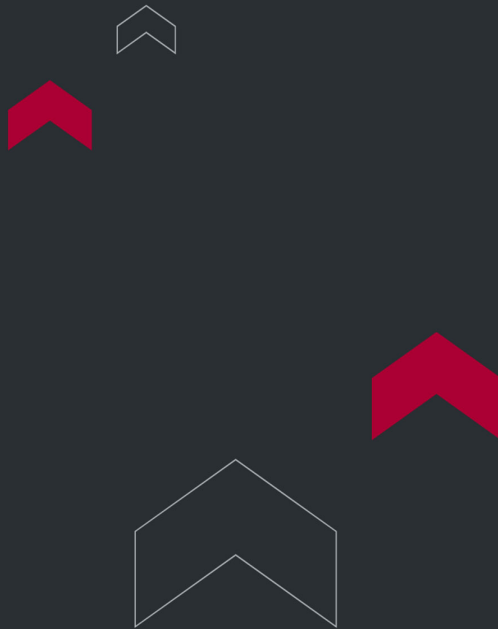

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

