



GOLD PURITY DENSITY TESTER TTD-214K

Gold Purity Density Tester TTD-214K specifically designed for precious metals, ideal for small and lightweight samples such as gold sand (dimensions: 148 x 100 x 85 mm).

Maximum Weight Capacity: 210g

Purity Precision: 100.0% or 100.00%

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

APPLICATIONS

Suitable for

Suitable for the jewelry industry, banks, pawnshops, and gold and precious metal research laboratories.

Principle

Using the buoyancy method based on the Archimedean principle, gold karat, density, and purity percentage can be displayed quickly and accurately.

Mode

1. Gold Purity Measuring Mode - Measures gold purity.
2. Platinum Purity Measuring Mode - Measures platinum purity.
3. Alloy Purity and Density Measuring Mode 1 - Measures the purity of other precious metals.
4. Alloy Purity and Density Measuring Mode 2 - Measures the purity of other precious metals.

FEATURES


- Gold Purity Density Tester TTD-214K equipped with RS-232 and USB interfaces.
- Infrared temperature sensors automatically detect and compensate for water temperature.
- Provides quick and accurate determination of precious metal purity percentage and gold karat number in just two measurement steps.
- Leaves no trace of contamination after measurement, eliminating the need for sulfuric acid solution or standard test samples.
- Special software design directly displays gold karat number and purity percentage (including gold-copper-silver, gold-silver, and gold-copper alloys).
- Features a setting function for the lower limit of karat value, allowing the karat number to be freely set from 1 to 24K.
- Includes a function for setting the lower limit of PT value, with PT value adjustable from 1 to 1000.
- Capable of measuring the purity of other precious metals, such as platinum, silver, palladium, copper, rhodium, etc.
- Unable to test samples with gems and hollows. Can be utilized as a general balance.

- Utilizes a large tank design to reduce errors caused by the buoyancy of the hanging rail. The water tank size is 148 × 100 × 85mm.

TECHNICAL SPECIFICATIONS

Model		TTD-214K
Max Weight		210g
Weighing Precision		0.0001g
Density Precision		0.0001g/cm ³
Gold Mode	Karat	0.1K and 0.01K display (switchable)
	Purity	0.1% and 0.01% display (switchable)
PT Mode	PT	1PT and 0.1PT display (switchable)
	Purity	0.11% and 0.01% display (switchable)
Alloy Mode		0.1% and 0.01% display (switchable)

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

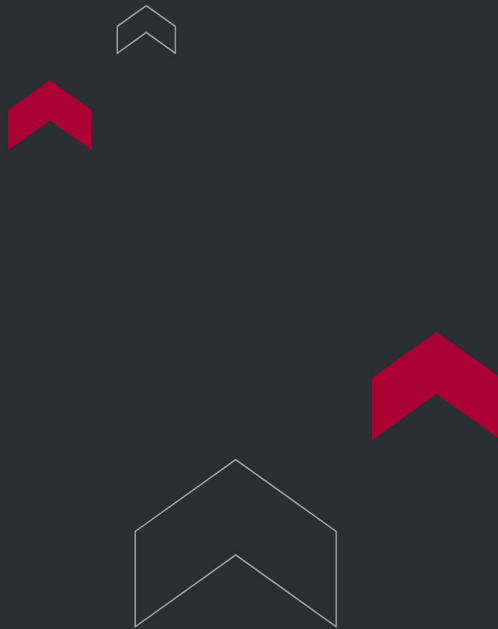

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 🇨🇦

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 🇺🇸

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 🇦🇪

Phone: +971 488 19252
Email: me@torontech.com

