



SOLID MATERIALS DENSITY TESTER TTD- 300E/600E

Solid Materials Density Tester TTD-300E/600E specifically designed for solid rubber materials, suitable for general-size samples with dimensions of 148mm in length, 100mm in width, and 85mm in height.

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Maximum Weight Capacity: Supports a maximum weight of either 300g or 600g.

Precision in Density Measurement: Offers density measurement precision up to 0.001g/cm³.

APPLICATIONS

Suitable for

Laboratory of solder, brass, bronze, silver, nickel, zinc, plastic, rubber, wire and cable, glass, film, hard alloy, soft synthetic leather, leather, elastic material, oil seal, oil ring, new materials, and other recycling industries.

Principle


Referring to the standards of ASTM D297, D792, JIS K6530, ISO2781, 1183and by applying the buoyancy of the object immersed in the liquid and hydrostatics principle, the density, volume, and mixture ratio of dense solid materials can be shown.

FEATURES

- Density Measuring Mode (DS): Measures both density and volume.
- Mix Ratio Density Measuring Mode (Mix1): Determines the proportion of the primary material.
- Mix Ratio Density Measuring Mode (Mix2): Determines the proportion of the secondary material.

TECHNICAL SPECIFICATIONS

Model	TTD-300E	TTD-600E
Max Weight	300g	600g
Weighing Precision	0.01g	
Density Precision	0.001g/cm3	
Volume Precision	0.001g/cm3	
Mixing Ratio 1	0.01%	
Mixing Ratio 2	0.01%	

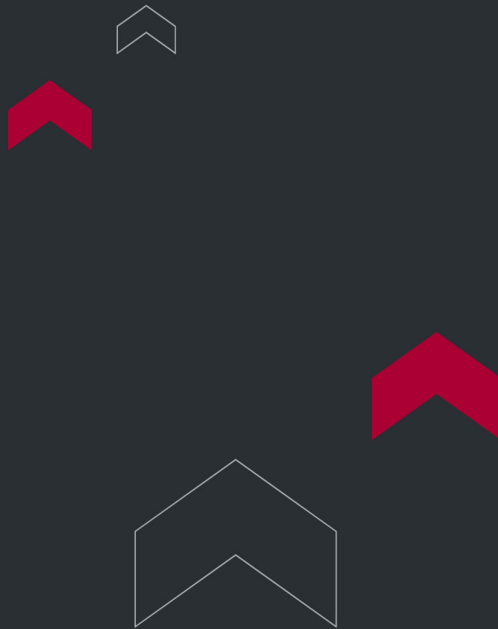

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

