



LPG COPPER STRIP CORROSION

LPG Copper Strip Corrosion TT-1838

LPG Copper Strip Corrosion TT-1838 tester adheres to ASTM D1838, the Standard Test Method for Copper Strip Corrosion by Liquefied Petroleum (LP) Gases.

Copper corrosion limits ensure that copper and copper-alloy fittings and connections, commonly utilized in various types of equipment for utilization, storage, and transportation, do not deteriorate. This apparatus is designed to detect the presence of components in liquefied petroleum gases that may cause corrosion to copper.

APPLICATIONS

The primary application of this tester is to detect the presence of components within Liquefied Petroleum (LP) Gases that are corrosive to copper.

This is done to:


- Ensure that the LP gas meets quality standards as defined by the ASTM D1838 test method.
- Prevent the deterioration of copper and copper-alloy fittings and connections that are used in equipment for storing, transporting, and utilizing LP gases.

FEATURES

- Utilizing a digital display temperature controller for highly accurate temperature control.
- The apparatus includes two pairs of experimental bombs, enhancing experimental efficiency.
- Equipped with a corrosion-resistant stainless steel bath.
- The appearance is designed with decency and beauty in mind.

TECHNICAL SPECIFICATIONS

1. Rated Voltage: AC 220V $\pm$ 10% 50Hz or 110V $\pm$ 10% 60Hz
2. Power: 2000W
3. Temperature Control: Digital display temperature controller
4. Temperature Range: Room temperature ~ 90°C
5. Control Precision: $\pm$ 0.5°C
6. Temperature Sensor: Pt100 (Platinum resistor)
7. Bath Holes: Two
8. Mixing Method: Motor stirring
9. Ambient Requirements: Temperature 10~40°C; Humidity $\leq$ 85%

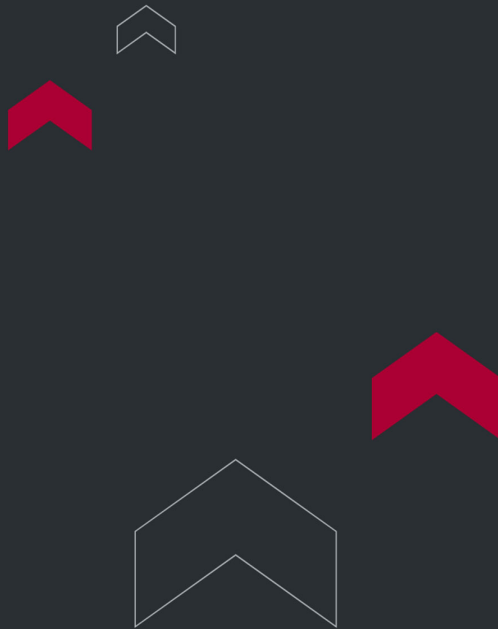

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 🇨🇦

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 🇺🇸

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 🇦🇪

Phone: +971 488 19252
Email: me@torontech.com

