

DENSITY TEST BY API HYDROMETER



Density Test by API Hydrometer TT-1298

Density Test by API Hydrometer TT-1298 conforms to ASTM D1298 Standard Test Method for Density, Relative Density, or API Gravity of Crude Petroleum and Liquid Petroleum Products by Hydrometer Method.

Density Test by API Hydrometer TT-1298 facilitates laboratory determination using a glass hydrometer along with a series of calculations to determine the density, relative density, or API gravity of crude petroleum, petroleum products, or mixtures of petroleum and nonpetroleum products typically handled as liquids. These measurements are made at existing temperatures and then corrected to 15°C or 60°F through a series of calculations and international standard tables.

APPLICATIONS

The primary application of this apparatus is for the laboratory determination of density, relative density, or API gravity of various liquid products, in compliance with the ASTM D1298 standard.

It is specifically used for testing:

- Crude petroleum
- Petroleum products
- Mixtures of petroleum and non-petroleum products that are handled as liquids




FEATURES

- The apparatus can conduct two-way experiments, significantly enhancing test efficiency.
- Featuring a stainless steel heater, the apparatus is corrosion-resistant, ensuring a long working life.
- A digital display meter is used to show temperature, providing high precision temperature control.
- Simple operation is facilitated by an elegant structure.

TECHNICAL SPECIFICATIONS

1. Rated voltage: AC 220V $\pm 10\%$ 50Hz (110V also available)
2. Power: 2000W
3. Bath temperature: 20°C to 95°C
4. Graduate flask: 1000ml
5. Temperature control: Digital display temperature controller
6. Temperature sensor: Pt100
7. Precision of temperature control: $\pm 0.2^\circ\text{C}$
8. Slot holes: Two holes
9. Mixing method: Mixing motor, 1200r/min
10. Ambient requirement: Temperature 10°C to 40°C; Humidity $\leq 85\%$

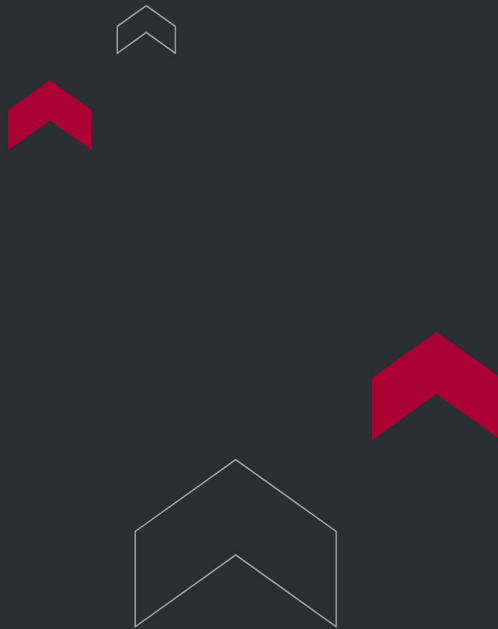

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

