



DIGITAL CONDUCTIVITY METER

Digital Conductivity Meter TT-2624

Conductivity Tester for Aviation and Distillate Fuels TT-2624 is designed according to the ASTM D2624 Standard Test Methods for Electrical Conductivity of Aviation and Distillate Fuels.

These methods encompass the determination of electrical conductivity in aviation and distillate fuels, both with and without a static dissipator additive. Typically, these test methods provide a measurement of conductivity when the fuel is uncharged, or electrically at rest, commonly referred to as the rest conductivity.

APPLICATIONS

The primary application of this tester is to determine the electrical conductivity (rest conductivity) of aviation and distillate fuels, in compliance with the ASTM D2624 standard.

It is specifically used for:

- Aviation fuels (with or without static dissipator additives).
- Distillate fuels (with or without static dissipator additives).
- Other light oil products.
- Environmental protection fuel additive detection.

FEATURES

- Technology: This tester uses single-chip microcomputer photoelectric isolation technology for high stability and repeatability.
- Key Features: It has a short polarization time, extended sampling time, and high accuracy.

- **Functionality:** The device includes an integrated temperature measurement function.
- **Application:** It's primarily used for measuring the conductivity of light oil products, but it can also be used for environmental protection fuel additives.
- **Included Components:** The configuration comes with a temperature measuring probe, a conductivity pool, and a detector.
- **Measurement Range:** The conductivity measurement range is from 1 to 2000 $\mu\text{S/m}$ with a resolution of 1 $\mu\text{S/m}$.


TECHNICAL SPECIFICATIONS

1. Conductivity Test:

- 1) Range: 0 to 1999 $\mu\text{S/m}$
- 2) Resolution: 1 $\mu\text{S/m}$
- 3) Accuracy: $< 2\% \pm 1 \mu\text{S/m}$
- 4) Measurement error: Less than the requirement of ASTM D2624
- 5) Ambient temperature: 0°C to 40°C
- 6) Relative humidity: $\leq 85\% \text{ RH}$
- 7) Power: 9V
- 8) Dimensions: 80 × 140 × 35 mm
- 9) Calibration value: Certificate value $\pm 2 \mu\text{S/m}$

2. Temperature Test:

- 1) Temperature range: -50°C to 50°C
- 2) Resolution: 0.1°C
- 3) Accuracy: $< 1^\circ\text{C}$

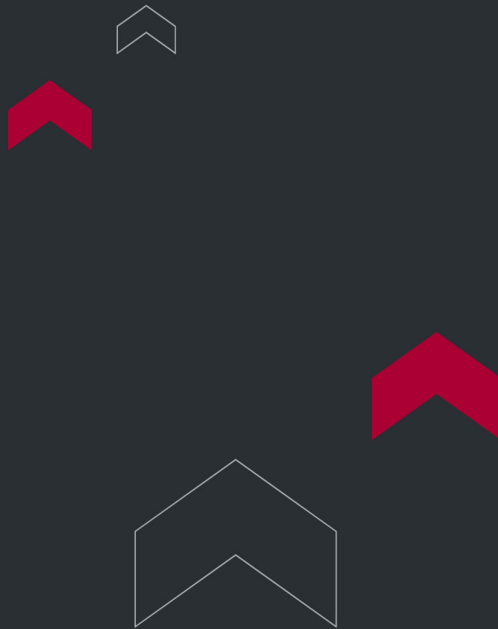

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

