



TORONTECH



FULLY AUTOMATED COLLECTION CUP WASHER

Fully Automated Collection Cup Washer is a cost-effective and user-friendly device designed to gently mix drugs in collection cups, dissolving active ingredients in solvents for analysis.

APPLICATIONS

- USP <601> Aerosols, Nasal Sprays, Metered-Dose Inhalers, and Dry Powder Inhalers
- European Pharmacopoeia (Ph. Eur.) 2.9.18: Preparations for Inhalation – Aerodynamic Assessment of Fine Particles
- Chinese Pharmacopoeia (ChP) 2020 – Inhalation Product Testing

Fully Automated Collection Cup Washer


FEATURES

1. Compatible with multi-stage cascade impactors, offering convenience and cost-efficiency.
2. Ensures consistent, accurate, and repeatable drug dissolution.
3. Features smart control, customizable settings, simple operation, and a clear display.

TECHNICAL SPECIFICATIONS

Dimensions	800x200x180mm
Weight	7.3KG
Voltage	220V, 50Hz (110V is also available)

Power Rate	120W
Max. Angle	20°
Max. Speed	20°/s
Running Accuracy	±1°
Max. Runs	30000

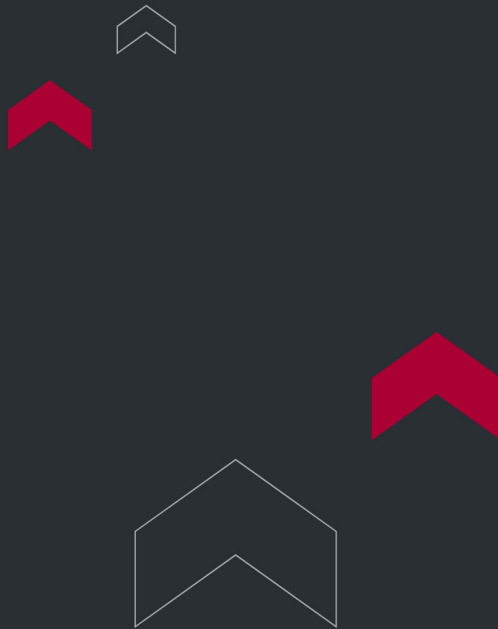

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

