



TORONTECH



SURGICAL TABLE - FULLY ELECTRIC MOBILE TORONSURGE™ -1000B

Surgical Table – Fully Electric Mobile ToronSurge™-1000B mobile operating table is designed for maximum versatility and ease of use across a wide range of surgical procedures. Built with high-quality 304 stainless steel, the bed and column are durable, easy to clean, and resistant to contaminants—ensuring superior hygiene and longevity.

APPLICATIONS

The table surface is divided into five sections: headrest, backrest, seat, and two detachable leg sections, offering exceptional flexibility for various surgical positions.

With 340mm horizontal translation, this table provides an optimal imaging area for C-arm procedures and is compatible with X-ray film boxes, making it an excellent choice for modern operating rooms.

This multi-functional electric surgical table is ideal for general surgery, abdominal procedures, obstetrics, gynecology, ENT, urology, anorectal surgery, and orthopedics.

FEATURES

Key Features

1. Double-Joint Headrest for Precision Positioning

The adjustable double-joint headrest allows precise height and angle adjustments, ensuring patient comfort and optimal surgical access.

2. Electrically or Gas Spring-Controlled Leg Plates

The leg plates can be manually adjusted with a gas spring or upgraded to full electric control for enhanced convenience and efficiency.

3. Optional Dual-Control System

The handheld remote control and optional panel control provide redundant safety measures, ensuring continuous operation in any scenario.

4. Heavy-Duty German TENTE Casters

High-quality wear-resistant and silent casters offer smooth mobility while maintaining maximum stability.

5. Built-In Rechargeable Battery

The high-performance rechargeable battery supports up to 50 surgeries on a single charge, with an AC power backup for uninterrupted operation.

6. Advanced Customization Options

Optional one-button reset, positive flexion, and reverse flexion functions enhance surgical adaptability and ease of use.

TECHNICAL SPECIFICATIONS

Model	Surgical Table - Fully Electric Mobile ToronSurge™-1000B
Dimensions (L × W)	2100mm × 550mm
Height Adjustment	1000mm - 700mm
Head Plate Adjustment	+45° / -90°
Back Plate Adjustment	+90° / -17°
Leg Plate Adjustment	+15° / -90° / -90°

Trendelenburg / Reverse Trendelenburg	25° / 25°
Lateral Tilt (Left & Right)	15° / 15°
Kidney Bridge Elevation	≥110mm
Horizontal Sliding	345mm
Flex / Reflex Function	Combination Operation
X-ray Compatibility	Optional X-ray Board
Control System	Standard Panel & Optional Hand Remote
Motor System	Electric Motor
Voltage Options	220V / 110V
Frequency	50Hz / 60Hz
Power Capacity	1.0 KW
Battery Backup	Yes
Mattress Type	Memory Foam Mattress
Material	304 Stainless Steel
Maximum Load Capacity	200 KG
Warranty	1 Year

Standard Accessories

Accessories for Enhanced Functionality



Anaesthesia Screen

Provides a privacy barrier during procedures.



Fixing Clamp

Securely attaches accessories to the table rail.



Stomach Support Rod

Supports abdominal positioning during surgery.




IV Stand

Holds IV bags for fluid administration.



Fixing Clamp

Securely attaches accessories to the table rail.

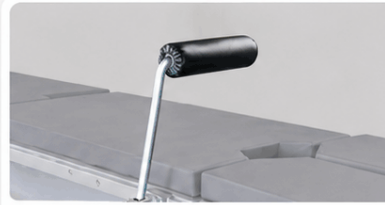

Knee Support Rod

Provides support and stability for the knee.




Arm Rest (PU)

Offers comfortable support for the arm.




Fixing Clamp

Securely attaches accessories to the table rail.

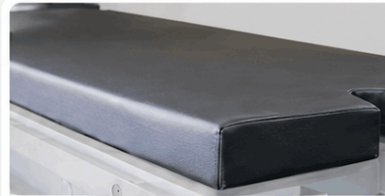

Seamless Memory Pad

Conforms to the body for enhanced comfort.

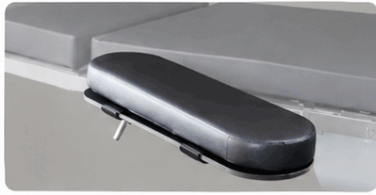


Arm Rest (Padded)

Offers comfortable support for the arm.



Fixing Clamp

Securely attaches accessories to the table rail.

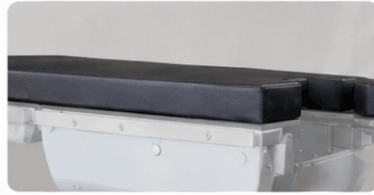

Memory Pad

Conforms to the body for enhanced comfort.




Arm Rest (Padded)

Offers comfortable support for the arm.




Fixing Clamp

Securely attaches accessories to the table rail.

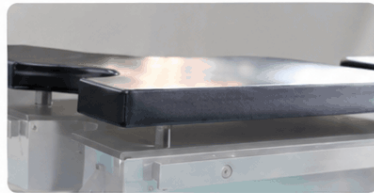

PU Pad

Durable surface for added support.

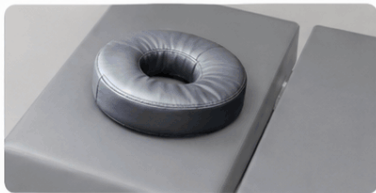

Head Ring

Provides head support and stability.



Body Belt (Buckle Type)

Secures the body for safety and stability.



Waste Basin

For convenient collection of waste.



Lateral Support Pad

Provides lateral support and positioning stability.



Body Belt

Secures the body for safety and stability.



Operation Support

Enhances support and positioning during procedures.



Prone Position Pad

Supports the body in prone position for procedures.



Leg Belt

Secures the legs to help maintain proper positioning.




Touchscreen Hand Remote

Intuitive touchscreen control for table adjustments.



Stomach Support Pad

Supports the abdomen to improve patient comfort.

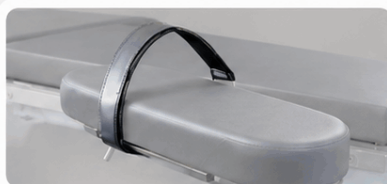

Arm Belt

Secures the arms to reduce movement and enhance safety.



Button Hand Remote


Easy-to-use button control for table positioning.

Accessories for Surgical Table

**Leg Support
(Padded)**



Wrist Belt

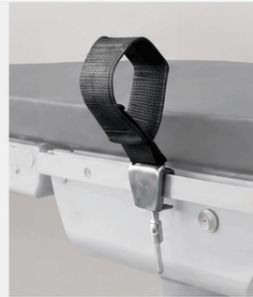

**Hand
Remote**



**Leg Support
(PU)**



Foot Plate



**Lateral
Arm Support**



**External
Kidney Bridge**



**Shoulder
Support**



**Instrument
Tray**



Arm Surgery
Stand



Body
Support



Orthopedic
Traction 10 05



Clamp



Body
Support




Orthopedic
Traction 10 06



Clamp




Head Rest

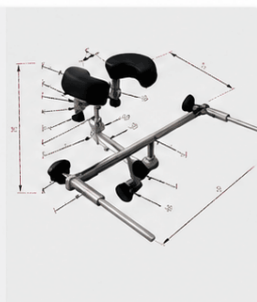

Three Point
Skull Clamp



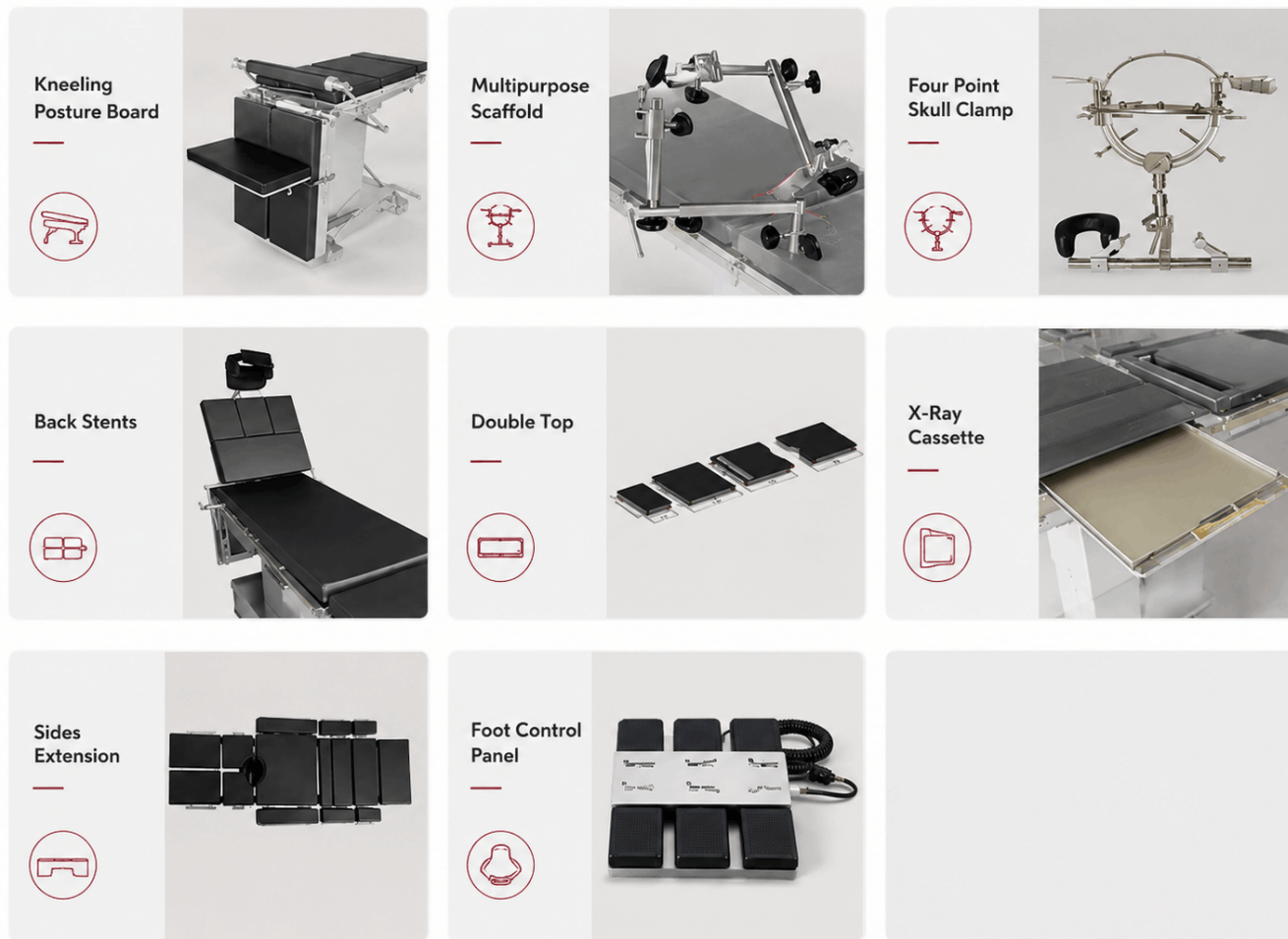

Accessory	Quantity
Anesthesia Screen	1 piece
Body Support	1 pair
Arm Supports	1 pair
Shoulder Supports	1 pair
Leg Supports	1 pair

Accessory	Quantity
Foot Supports	1 pair
Kidney Bridge Handle	1 piece
Memory Foam Mattress	1 set
Fixing Clamps	8 pieces
Long Fixing Clamps	1 pair
Hand Remote	1 piece
Power Cable	1 piece



Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

