



TORONTECH



SURGICAL TABLE - HYDRAULIC TORONSURGE™ -3000Y

ToronSurge™ -3000Y hydraulic surgical table offers precise manual adjustments with pneumatic and hydraulic support for seamless operation. Designed for versatility, it accommodates a range of procedures, ensuring stability, mobility, and surgeon comfort with its ergonomic T-shaped base. Its durable, antistatic foam mattress enhances hygiene and patient safety.

APPLICATIONS

- **General Surgery** - Provides stable manual positioning with hydraulic and pneumatic support for chest and abdominal procedures.
- **Gynecology & Obstetrics** - Adjustable positioning and ergonomic design to support gynecological examinations and childbirth-related surgeries.
- **Urology & Anorectal Surgery** - Facilitates precise patient positioning and stability for lower abdominal and anorectal treatments.
- **Orthopedic Surgery** - Supports secure limb positioning and multiple angle adjustments for fracture management and orthopedic interventions.
- **Emergency & Trauma Surgery** - Manual hydraulic adjustments ensure reliable operation even without external power sources.
- **Imaging & Diagnostics** - Radiolucent design and optional 360° rotation enable efficient use with X-ray and C-arm equipment.

FEATURES

Key Features of Surgical Table - Hydraulic ToronSurge™ -3000Y

- 1. All positions are manually adjustable and supported by pneumatic and hydraulic systems, ensuring effortless operation, precision, and reliability.
- 2. This operating table supports various procedures, including chest and abdominal surgeries, gynecology, urology, anorectal treatments, and orthopedics.
- 3. The seamless, high-quality foam mattress is antistatic and resistant to disinfectants, enhancing durability and hygiene.
- 4. Engineered for both stability and mobility, the T-shaped base features an ergonomic design that provides ample foot space for surgeons.

TECHNICAL SPECIFICATIONS

Model	ToronSurge™ -3000Y
Length/Width	2010/500mm
Lowest/Highest height	760/1010mm
Trendelenburg/Reverse Trendelenburg	≥22°/22°
Lateral Tilt (Left/Right)	≥20°/20°
Head Plate Up/Down	≥60°/70°
Leg Plate Down/Outward/Up	≥15°/90°/85°
Back Plate Up/Down	≥75°/15°
Kidney Bridge Elevation	>120mm Optional

Model	ToronSurge™ -3000Y
Leg Plate adjustment	Manual
Double Deck	Yes
360° Rotation	Optional

Standard Accessories

Accessories for Enhanced Functionality



Anaesthesia Screen

Provides a privacy barrier during procedures.



Fixing Clamp

Securely attaches accessories to the table rail.



Stomach Support Rod

Supports abdominal positioning during surgery.




IV Stand

Holds IV bags for fluid administration.



Fixing Clamp

Securely attaches accessories to the table rail.

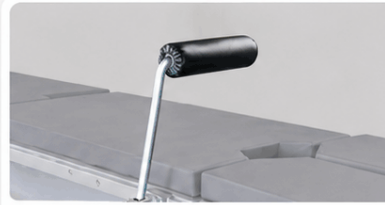

Knee Support Rod

Provides support and stability for the knee.




Arm Rest (PU)

Offers comfortable support for the arm.




Fixing Clamp

Securely attaches accessories to the table rail.

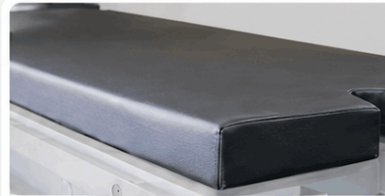

Seamless Memory Pad

Conforms to the body for enhanced comfort.

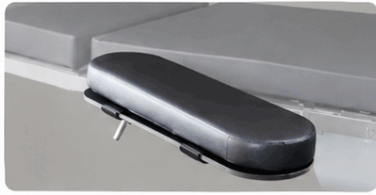


Arm Rest (Padded)

Offers comfortable support for the arm.



Fixing Clamp

Securely attaches accessories to the table rail.

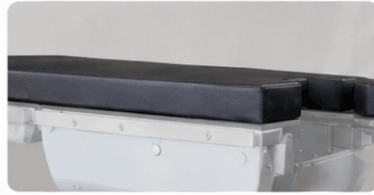

Memory Pad

Conforms to the body for enhanced comfort.




Arm Rest (Padded)

Offers comfortable support for the arm.




Fixing Clamp

Securely attaches accessories to the table rail.

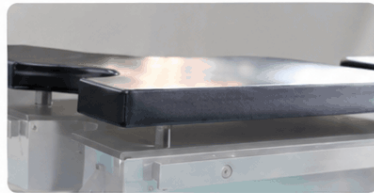

PU Pad

Durable surface for added support.

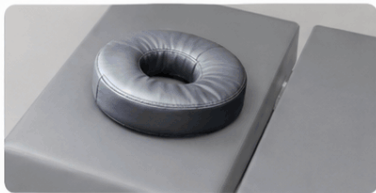

Head Ring

Provides head support and stability.



Body Belt (Buckle Type)

Secures the body for safety and stability.



Waste Basin

For convenient collection of waste.



Lateral Support Pad

Provides lateral support and positioning stability.



Body Belt

Secures the body for safety and stability.



Operation Support

Enhances support and positioning during procedures.



Prone Position Pad

Supports the body in prone position for procedures.



Leg Belt

Secures the legs to help maintain proper positioning.




Touchscreen Hand Remote

Intuitive touchscreen control for table adjustments.



Stomach Support Pad

Supports the abdomen to improve patient comfort.

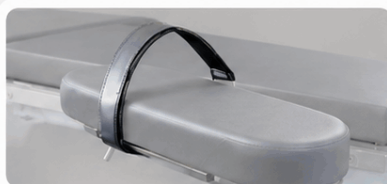

Arm Belt

Secures the arms to reduce movement and enhance safety.



Button Hand Remote


Easy-to-use button control for table positioning.

Accessories for Surgical Table

**Leg Support
(Padded)**



Wrist Belt

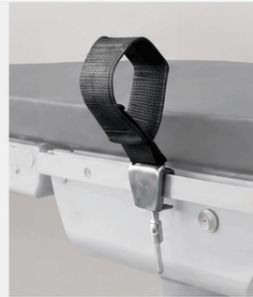

**Hand
Remote**



**Leg Support
(PU)**



Foot Plate



**Lateral
Arm Support**



**External
Kidney Bridge**



**Shoulder
Support**



**Instrument
Tray**



Arm Surgery
Stand



Body
Support



Orthopedic
Traction 10 05



Clamp



Body
Support




Orthopedic
Traction 10 06



Clamp




Head Rest

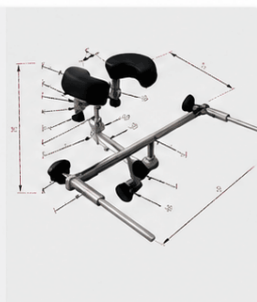

Three Point
Skull Clamp



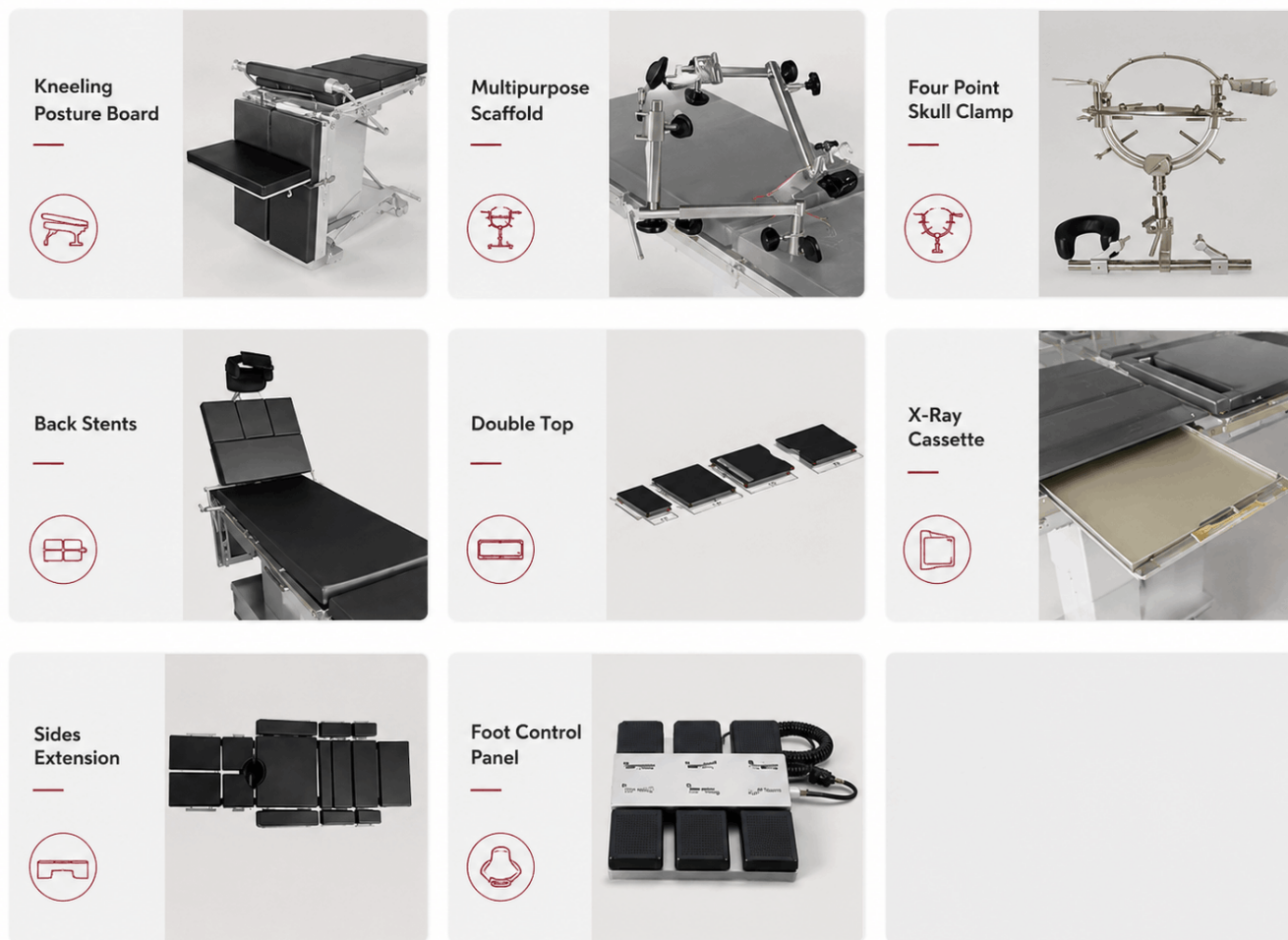

Accessory	Quantity
Anesthesia Screen	1 piece
Arm Support	1 pair
Body Support	1 pair
Leg Support	1 pair
Mattress	1 set



Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

