



TORONTECH

MEDICATION CART (STANDARD)



Medication Cart (Standard) - ToronCAB™ 1230

ToronCAB™ 1230 standard medication cart is a versatile medical cart that organizes medicines, supplies, and sharps disposal for bedside treatment. Two lockable drawers, secured by a single central lock, hold medications and consumables, while a medicine cassette holder keeps dose trays sorted and accessible. A sliding work surface gives staff room to prepare and administer treatments at the point of care. Four braking swivel casters move the cart between rooms and lock it steady during use.

Nurses and clinical staff use the Torontech ToronCAB™ 1230 for medication rounds and routine treatments across hospital wards and clinics. A dedicated needle disposal container and two waste bins let staff drop sharps and general waste into the right container at the moment of use, supporting safe segregation. A disinfectant holder keeps hand or surface disinfectant within reach for hygiene between patients. The removable PVC top mat and molded surfaces wipe down quickly, which supports infection-control routines on busy units.

For North American facilities, the central-locking drawers help keep medications secured between rounds, supporting controlled-access handling of drugs at the unit level. The compact L620 by W500 mm footprint moves easily through doorways and fits in tight treatment rooms and beside beds. Designed for daily clinical use, the ToronCAB™ 1230 brings medicine storage, a work surface, sharps disposal, and disinfectant together on one mobile platform, so staff handle more of each task at the bedside.

APPLICATIONS

Medication Cart (Standard) - ToronCAB™ 1230 Applications

The ToronCAB™ 1230 supports medication delivery, routine treatment, and safe sharps handling across departments that administer drugs at the bedside. Its drawers, medicine holder, and disposal containers fit a range of clinical settings.

- Hospital Inpatient Wards: carries medicines and supplies for scheduled drug rounds, so nurses dispense at the bedside without returning to the medication room.
- Treatment and Procedure Rooms: provides locked drug storage, a work surface, and sharps disposal for injections, dressing changes, and minor procedures.
- Outpatient Clinics and Day Units: moves medicines and consumables between treatment bays for back-to-back patient appointments.
- Intensive Care and High-Dependency Units: keeps frequently used drugs, disposables, and sharps containers organized and mobile in high-activity areas.
- Vaccination and Injection Clinics: organizes dose trays in the medicine cassette holder and disposes of needles immediately in the on-board container.
- Long-Term Care and Rehabilitation Facilities: supports medication rounds and treatments with secure storage and disinfectant on hand in shared environments.
- Recovery and Observation Areas: positions medicines, a work surface, and waste handling near staff for fast care after procedures.



FEATURES

Medication Cart (Standard) - ToronCAB™ 1230 Key Features

The ToronCAB™ 1230 combines secure medicine storage, organized dose handling, and safe waste disposal in a compact medication cart for daily clinical use.

- Two Lockable Drawers: drawers store drugs and consumables and stay organized at working height for quick access during rounds.
- Single Central Lock: one lock secures both drawers at once, so staff close down medications in a single action when leaving the cart.
- Medicine Cassette Holder: a dedicated holder arranges dose trays in view, so staff select the right medication without searching a drawer.
- Sliding Work Surface: a pull-out table extends for medication preparation and retracts to keep the cart compact in transit.
- Needle Disposal Container: an on-board sharps container lets staff dispose of used needles immediately and safely at the bedside.
- Dual Waste Bins: two side-mounted bins separate general and clinical waste at the point of use.
- Disinfectant Holder: keeps disinfectant within reach so staff clean hands and surfaces between patients.
- Four Individual-Brake Casters: each swivel caster locks on its own for a stable base during use and free movement between rooms.
- Corner Anti-Bumper Guards: bumpers on all four corners cushion contact with walls and doorframes, protecting the facility and the cart.
- Removable PVC Top Mat: the protective mat lifts off for cleaning and shields the work surface from spills and scratches.

THEORY & METHOD

Theory and Method

The ToronCAB™ 1230 organizes medication delivery around secure drawers and an open medicine cassette holder. Two drawers hold drugs and consumables behind a single central lock, so staff secure everything in one action when leaving the cart. The medicine cassette holder arranges dose trays or cassettes in view, so staff find the right medication quickly instead of searching a deep drawer. Working at this height keeps frequently used items reachable without bending.

Safe disposal and hygiene shape the rest of the cart's design. A needle disposal container sits on the cart so staff drop used sharps straight into a sealed container at the bedside, rather than carrying them to a distant bin. Two waste bins separate general and clinical waste at the point of use, which keeps segregation simple during procedures. A disinfectant holder keeps cleaning solution at hand, so staff disinfect surfaces and hands between patients without stepping away from the cart.

Mobility comes from four swivel casters, each with its own brake. Locking all four casters holds the cart steady while staff prepare or administer medications, and releasing them gives free movement in any direction through wards and doorways. The sliding work surface extends when staff need room to work and retracts to keep the compact cart easy to maneuver. Corner anti-bumpers absorb knocks during transport, so the cart and the building avoid damage over years of daily use.

TECHNICAL SPECIFICATIONS


Medication Cart (Standard) - ToronCAB™ 1230 Technical Specifications

Parameter	Specification
External Dimensions	L620 × W500 × H1020 mm (±10 mm)
Medication Storage	Two lockable drawers with a single central lock
Medicine Organization	Medicine cassette holder for dose trays

Parameter	Specification
Work Surface	Molded top with removable PVC mat, plus a sliding side table
Waste and Sharps	Two waste bins plus a needle disposal container
Hygiene	Disinfectant holder
Mobility	Four swivel casters with individual brakes

Standard Configuration

Item	Quantity	Description
Drawer	2 pcs	Lockable drawers for medications and consumables
Central Lock	1 pc	Single lock securing both drawers in one action
Medicine Cassette Holder	1 pc	Holder for arranging dose trays and cassettes
Sliding Working Table	1 pc	Pull-out work surface for medication preparation
Needle Disposal Container	1 pc	On-board sharps container for used needles
Waste Bin	2 pcs	Side-mounted bins for general and clinical waste
Disinfectant Holder	1 pc	Holder keeping disinfectant within reach
Top Mat (PVC)	1 pc	Removable PVC mat protecting the top work surface
Anti-Bumper	4 pcs	Corner-mounted guards protecting walls, fixtures, and the cart
Individual Brake Caster	4 pcs	Swivel casters, each with a dedicated brake for stable positioning

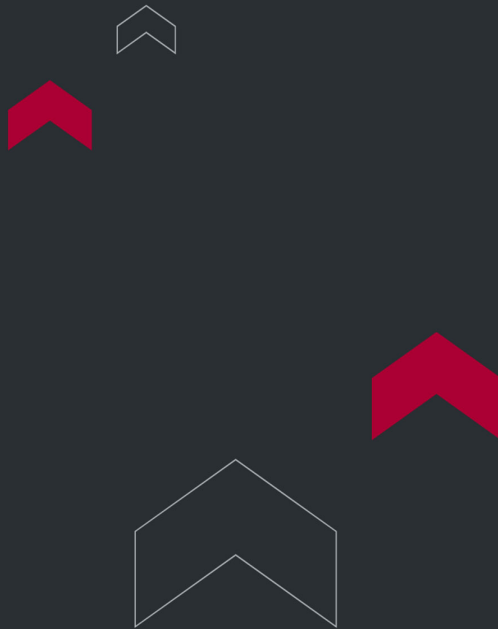

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

