



CRASH CART (PREMIUM)

Crash Cart (Premium) - ToronCAB™ 1460

ToronCAB™ 1460 is a premium medical crash trolley that brings organized drug storage and resuscitation equipment to the bedside during a code. Four drawers hold emergency medications across sixteen adjustable compartments, and a lockable tabletop medicine box with a light-proof cover protects sensitive drugs on the work surface. An adjustable defibrillator tray, CPR board, and oxygen cylinder holder keep resuscitation equipment with the cart, while a height-adjustable IV pole supports infusions. Four braking 125 mm casters move the cart and lock it steady.

Hospital teams use the Torontech ToronCAB™ 1460 in settings that demand both fast response and a durable, well-organized cart. An impact-resistant HDPE body resists knocks and cleans easily, while a welded steel internal frame keeps the structure rigid under daily use. Four drawers run on three-step, self-locking silent slides, and a lateral storage box adds room for bulkier supplies. An ergonomic HDPE push handle and 125 mm silent casters let one person move the loaded cart smoothly and quietly through corridors.

For North American facilities, a central control lock and disposable breakaway seals secure all drawers in one action and signal at a glance whether the cart is sealed and fully stocked between codes. Manufactured under ISO 13485 and ISO 9001 quality systems, the cart suits departments that want a certified, premium emergency platform. Designed for code response, the ToronCAB™ 1460 brings classified drug storage, a defibrillator tray, oxygen, a CPR board, and power together, so the team works from one organized, ready station.

APPLICATIONS

Crash Cart (Premium) - ToronCAB™ 1460 Applications

The ToronCAB™ 1460 supports emergency and resuscitation response in departments that want a premium, well-organized crash cart with certified construction. Its classified drug storage, defibrillator tray, oxygen, and CPR board fit demanding

acute settings.

- Emergency Departments: stands ready for cardiac arrests and trauma resuscitations, with classified drug storage and a defibrillator tray organized for immediate use.
- Intensive Care Units: keeps emergency drugs, oxygen, and a defibrillator within reach for unstable patients in high-acuity care.
- Operating Rooms and Recovery: supports response to intraoperative and post-anesthesia emergencies alongside the surgical team.
- Hospital Inpatient Wards: serves as a mobile code cart that staff bring to the bedside the moment a patient deteriorates.
- Catheterization and Procedure Labs: stays on standby for cardiac and interventional procedures where complications demand fast, well-stocked response.
- Pharmacy-Managed Medication Areas: holds light-sensitive and controlled drugs in the lockable tabletop box and sealed drawers for organized, secure storage.
- Maternity and Labor Wards: keeps emergency drugs and equipment ready for obstetric and neonatal resuscitation.



Standards

The ToronCAB™ 1460 is manufactured under internationally recognized quality and medical-device management standards, and the manufacturer's certification scope covers medical trolleys.

- ISO 13485:2016 - Medical Devices: Quality Management Systems, Requirements for Regulatory Purposes

FEATURES

Crash Cart (Premium) - ToronCAB™ 1460 Key Features

01

Superior Casters

High-grade silent casters with 125 mm diameter for smooth movement and anti-winding performance.



02

Mute Slide

Three-section silent drawer slides run smoothly and quietly, with a self-locking function.



03

Classified Storage


The trolley includes 4 drawers, each with classification bins, medicine cups, and card slots for organized medicine storage.



04

Built-in Medicine Box

Built-in medicine box on the countertop with a flip-top cover that can be locked to help keep medicines secure.

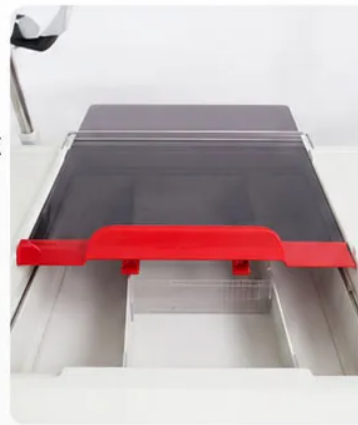

05

Durable Structure

The interior of the cart body uses a strong steel structure with firm welding for reliable stability.




06

Safety Bumper

HDPE blow-molded bumper construction provides toughness and helps protect the trolley from impact.



07

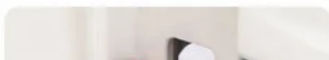

08



The ToronCAB™ 1460 combines premium construction, classified drug storage, and full resuscitation equipment in a hospital crash cart built for demanding code response.

- Lockable Tabletop Medicine Box: a built-in countertop box with a light-proof, locking cover protects light-sensitive and high-alert drugs on the work surface.
- Four Drawers with Sixteen Adjustable Compartments: drawers divide into resizable compartments, so facilities sort emergency drugs by class and adjust each space to fit.
- Adjustable Defibrillator Tray: a rotating tray positions the defibrillator at any angle and holds it with a fixed belt during transport.
- CPR Board: a light, high-strength backboard stores on the back and lifts out for firm chest compressions.
- Oxygen Cylinder Holder: a side holder carries an oxygen tank with the cart for immediate ventilation support.
- Height-Adjustable IV Pole: a telescoping pole with four hooks (rated five kilograms each) holds several infusions at once.
- Central Control Lock and Breakaway Seals: one lock secures all drawers, and disposable seals show at a glance whether the cart is sealed and fully stocked.
- Built-In Power Strip: on-board power runs the defibrillator or a monitor without depending on wall outlets.
- Impact-Resistant HDPE Body: a blow-molded body resists knocks and wipes down easily, with a welded steel internal frame for rigidity.
- Lateral Storage Box: a side box adds room for bulkier supplies beyond the drawers.
- Needle Disposal and Dual Waste Bins: an on-board sharps container and two waste bins handle sharps and waste safely at the point of use.
- Silent Self-Locking Slides and 125 mm Casters: three-step self-locking drawer slides and high-grade silent casters keep operation smooth and quiet.
- Ergonomic Push Handle: an HDPE push handle gives staff a firm, comfortable grip when moving a loaded cart.

THEORY & METHOD

Theory and Method

A crash cart performs only when staff can find each drug instantly, so the ToronCAB™ 1460 organizes medications for fast, classified access. Four drawers divide into sixteen adjustable compartments, so a facility sorts drugs by class and resizes each compartment to fit. A lockable tabletop medicine box keeps frequently grabbed or light-sensitive drugs on the work surface under a light-proof, locking cover. Organizing storage this way removes hesitation during a code while keeping medications secure and protected.

The cart carries the equipment a resuscitation demands, positioned for quick use. An adjustable defibrillator tray holds the defibrillator and rotates to face the team from any angle, with a fixed belt to keep the unit secure during transport. A CPR board on the back lifts out to place under the patient for firm compressions, and an oxygen cylinder holder carries a tank for ventilation. A height-adjustable IV pole adds four hooks rated for five kilograms each, so several infusions hang ready at once.

Build quality and security match the premium positioning. An HDPE blow-molded body gives the cart strong impact resistance and smooth, wipeable surfaces, while a welded steel internal frame holds everything rigid. A central control lock secures all drawers at once, and disposable breakaway seals lock the cart and stay intact until it is opened, giving a clear sign of readiness. Four 125 mm silent casters with brakes and an ergonomic push handle move the loaded cart quickly and lock it firmly at the bedside.

TECHNICAL SPECIFICATIONS

Crash Cart (Premium) - ToronCAB™ 1460 Technical Specifications

Parameter	Specification
External Dimensions	L965 × W560 × H980 mm (±10 mm)

Parameter	Specification
Body Construction	Blow-molded HDPE with welded steel internal frame
Medication Storage	Four drawers with sixteen adjustable compartments
Tabletop Medicine Box	Built-in box with lockable, light-proof cover
Tamper Evidence	Central control lock with disposable breakaway seals
Resuscitation Equipment	Adjustable defibrillator tray (with fixed belt), CPR board, oxygen cylinder holder
Infusion Support	Height-adjustable IV pole, four hooks (5 kg each)
Power	Integrated power strip
Additional Storage	Lateral storage box
Waste and Sharps	Two waste bins plus a needle disposal container
Drawer Slides	Three-step, self-locking silent slides
Casters	Four × 125 mm silent swivel casters with brakes
Handle	Ergonomic HDPE push handle

Standard Configuration

Item	Quantity	Description
Drawer	4 pcs	Drawers with sixteen adjustable compartments for classified drug storage
Central Lock	1 pc	Central control lock securing all drawers from one point
Disposal Plastic Seal	100 pcs	Disposable breakaway seals showing whether the cart has been opened
Built-In Tabletop Medicine Box	1 pc	Countertop box with a lockable, light-proof cover
I.V. Pole	1 pc	Height-adjustable pole with four hooks for infusions

Item	Quantity	Description
Defibrillator Shelf	1 pc	Adjustable tray holding the defibrillator, with a fixed belt
Needle Disposal Container	1 pc	On-board sharps container for used needles
Waste Bin	2 pcs	Side-mounted bins for general and clinical waste
Lateral Storage Box	1 pc	Side box for bulkier supplies
CPR Board	1 pc	Light, high-strength backboard for chest compressions
Oxygen Cylinder Holder	1 pc	Holder carrying an oxygen tank for ventilation
Power Strip	1 pc	Powers the defibrillator or a monitor at the cart
Individual Brake Caster	4 pcs	125 mm silent swivel casters, each with a brake

••• Complete Accessories •••

Thoughtful design for greater convenience and efficiency



CPR Backboard

CPR board on the back of trolley with light weight and high strength.



Power Strip & Oxygen Tank Holder

Power strip and oxygen tank holder on the side of trolley, convenient to use.



IV Pole with 4 Hooks

Height adjustable IV pole with four hooks, each hook's working load is 5 kgs.



Defibrillator Tray

Adjustable defibrillator tray, easy to place at any direction, with fixed belt.



Needle Disposal Container

Safe and convenient disposal of sharps to help prevent cross-infection.



Waste Bins (2 pcs)

Two separate waste bins for efficient waste classification and management.



Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

