



MEDICATION CART (ABS FEATURE)

Medication Cart (ABS Feature) - ToronCAB™ 1600

ToronCAB™ 1600 is an ABS medical treatment cart designed for dressing changes, wound care, and bedside procedures. Two railed ABS (Anti-lock Braking System) decks hold clean supplies and instruments, a drawer under the top deck stores smaller items, and a stainless steel basin sits ready for cleaning solutions or used materials. A waste bin under the second deck handles disposal at the point of care. Four 100 mm casters, two with brakes, move the cart between rooms and hold it steady during a procedure.

Nurses and clinical staff use the Torontech ToronCAB™ 1600 as a self-contained station for dressing and treatment rounds. Clean supplies stay in the drawer and on the decks, the basin holds antiseptic or collects soiled dressings, and the waste bin takes general waste, so a full procedure happens at one cart without back-and-forth. Molded ABS decks and a stainless steel frame wipe down quickly between patients, and corner bumpers protect walls and equipment during transport. The cart carries loads up to 100 kg.

The compact 720 by 460 mm footprint moves easily through corridors and into tight treatment areas, and the cart ships in knock-down form to lower freight cost. Built for daily use, the ToronCAB™ 1600 brings clean storage, a basin, and waste handling together for organized, mobile treatment.

APPLICATIONS

Medication Cart (ABS Feature) - ToronCAB™ 1600 Applications

The ToronCAB™ 1600 supports dressing, wound care, and minor treatment across departments that need clean storage, a basin, and waste handling on one cart. Its treatment-station layout fits a range of clinical settings.

- Treatment and Dressing Rooms: keeps clean supplies in the drawer, antiseptic in the basin, and waste in the bin for dressing changes and wound care.

- Hospital Inpatient Wards: brings a complete dressing station to the bedside during treatment rounds.
- Outpatient and Day Clinics: moves between treatment bays with everything a minor procedure needs in one place.
- Emergency Departments: supports wound cleaning and dressing during fast-moving care.
- Operating Room and Recovery: assists dressing and minor procedures alongside the surgical team.
- Dental Clinics: organizes tools, a basin, and waste for chairside procedures.
- Long-Term Care and Rehabilitation Facilities: supports routine wound care and dressing changes in shared environments.



FEATURES

Medication Cart (ABS Feature) - ToronCAB™ 1600 Key Features

The ToronCAB™ 1600 combines clean storage, a basin, and waste handling in a complete ABS treatment cart for dressing and wound care.

- Two Railed ABS Decks: two molded shelves with guardrails hold clean supplies and instruments and keep them from sliding off.
- Drawer Under the Top Deck: an ABS drawer stores smaller clean supplies within reach during a procedure.
- Stainless Steel Basin: a mounted basin holds antiseptic solution or collects soiled dressings, keeping the clean and dirty sides separate.
- Waste Bin Under the Second Deck: a plastic bin takes general waste at the point of care.
- Stainless Steel Frame with ABS Decks: smooth, wipeable surfaces stand up to cleaning between patients.
- 100 kg Load Capacity: the frame carries up to 100 kg across the decks for heavier tasks.
- Corner Bumpers: bumpers on each corner cushion contact with walls and equipment during transport.
- Four 100 mm Casters with Brakes: four casters roll smoothly over hospital floors, and two with brakes hold the cart steady when parked.
- Knock-Down Construction: the cart ships flat to lower freight cost and assembles on site without special tools.

THEORY & METHOD

Theory and Method

The ToronCAB™ 1600 organizes a dressing or treatment procedure around three zones on one cart. Clean supplies stay in the drawer and on the railed decks, within reach and protected from contamination. A stainless steel basin holds antiseptic solution or collects soiled dressings and instruments during the procedure, keeping the dirty side separate from the clean. A waste bin under the second deck takes general waste, so staff complete the clean, treat, and dispose steps without leaving the cart.


Construction supports hygiene and daily use. Molded ABS decks and a stainless steel frame give smooth, wipeable surfaces that stand up to cleaning between patients, and corner bumpers cushion knocks against walls and equipment. The frame carries loads up to 100 kg, and four 100 mm casters move the cart over hospital floors with brakes on two for stability. The cart assembles from a knock-down kit, which ships flat to cut freight cost and goes together on site.

TECHNICAL SPECIFICATIONS

Medication Cart (ABS Feature) - ToronCAB™ 1600 Technical Specifications

Parameter	Specification
External Dimensions	720 (L) × 460 (W) × 975 (H) mm (±10 mm)
Frame	Stainless steel
Decks	Two ABS engineering plastic decks with guardrails
Drawer	One ABS drawer under the top deck
Basin	Stainless steel basin
Waste Bin	Plastic waste bin under the second deck
Corner Protection	Bumpers on each corner
Casters	Four × 100 mm; two with brakes

Parameter	Specification
Load Capacity	100 kg
Construction	Knock-down (flat-pack)
Gross Weight	18 kg
Shipping Volume	0.15 m ³

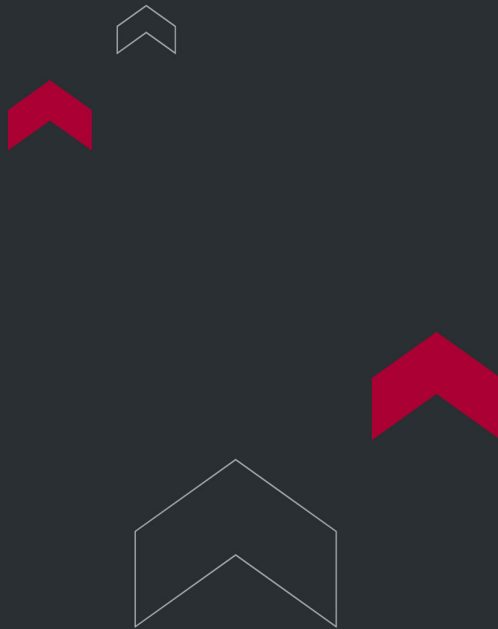

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 🇨🇦

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 🇺🇸

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 🇦🇪

Phone: +971 488 19252
Email: me@torontech.com

