



MEDICAL EXAM TABLES

Medical Exam Tables

Electric Exam Tables | IEC 60601-1 Safety-Aligned Positioning

Electric exam tables give clinics motorized height, backrest, and positioning control for exams and minor procedures. They are built around IEC 60601-1 safety practices for electric medical equipment. This range spans a base height-adjustable model, a mid-tier table with automated hygiene barrier, and a full multi-function table with tilt and translation. Physicians, nurses, and clinic staff use these tables across general practice and procedure rooms.

Each table shares a foot-pedal-controlled linear motor for height adjustment. Staff never need to bend or crank by hand. From there, the range splits by function. Some clinics need only reliable height adjustment. Others need added hygiene automation or full positioning control.

Why Choose Torontech for Your Medical Exam Tables

Torontech supplies testing and measurement equipment to laboratories and manufacturers across North America and worldwide, backed by responsive technical support and calibration service. The Medical Exam Tables range is built for reliable, standards-driven results, so you can standardize on one supplier as your testing needs grow. Our team helps you match the right configuration to your samples, applicable standards, and throughput, then supports you through installation and beyond.

Request a Quote

Tell us your patient volume, the positioning functions you rely on most, and your budget, and we will recommend the right table from the range. [Contact the Torontech team](#) for pricing and lead times.

Applications

Electric exam tables support a wide range of clinical settings where a patient needs fast, precise positioning.

- General Practice and Outpatient Clinics: foot-pedal height adjustment accommodates patients and clinicians of different sizes without manual cranking.
- Minor Procedure and Treatment Rooms: Trendelenburg tilt on the multi-function model supports procedures requiring a head-down position.
- Diagnostic Imaging Support: lateral tabletop translation positions a patient under fixed equipment without moving the whole table.
- High-Volume Outpatient Care: automated paper sheet changing reduces the time staff spend renewing the exam surface between patients.
- Multi-Provider Clinic Rooms: quick electric height changes reduce setup time between patients seen by different practitioners.

TECHNICAL SPECIFICATIONS

Medical Exam Table Range

The lineup spans three tiers, from a simple electric-lift table to a fully motorized multi-function platform.

Power Exam Table (ToronCare™ 1131)

The full multi-function model in the range. Four linear motors drive height, backrest, Trendelenburg and reverse Trendelenburg tilt, and a 300 mm lateral tabletop translation, plus an electric paper roll dispenser. It suits clinics that need imaging positioning support and procedure-specific tilt on one table.

Electric Examination Table (ToronCare™ 1132)


The base model in the range. A single linear motor handles foot-pedal height adjustment, and a gas spring drives the backrest through both an upward and downward range. It suits general practice rooms that need reliable, everyday height adjustment.

Electric Examination Table with Sheet Changer (ToronCare™ 1133)

A mid-tier model that adds a dedicated motor for automated disposable sheet changing, with full-width and half-width settings, on top of the base table's foot-pedal height adjustment and gas-spring backrest. It suits high-volume clinics that renew the exam surface often.

How to Choose the Right Medical Exam Table

- For imaging positioning or procedures needing tilt: choose the Power Exam Table (ToronCare™ 1131) for its lateral translation and Trendelenburg function.
- For straightforward, everyday exams on a budget: choose the Electric Examination Table (ToronCare™ 1132) for reliable height and backrest adjustment.
- For high patient turnover and hygiene priorities: choose the Electric Examination Table with Sheet Changer (ToronCare™ 1133) for automated barrier renewal.

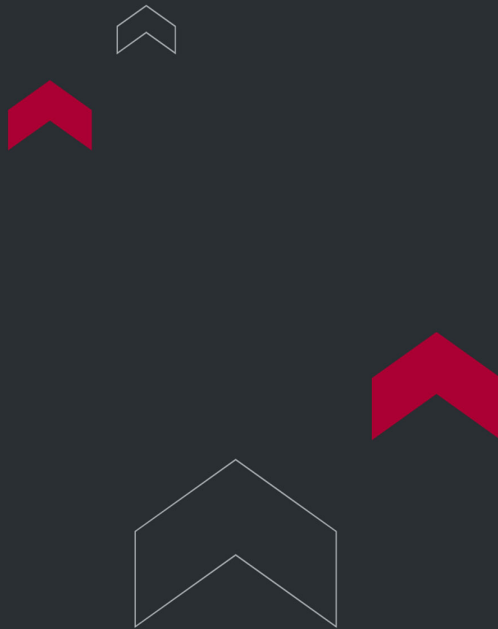

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

